




**MedLED<sup>®</sup>** Spectra  
Cordless Surgical Headlight

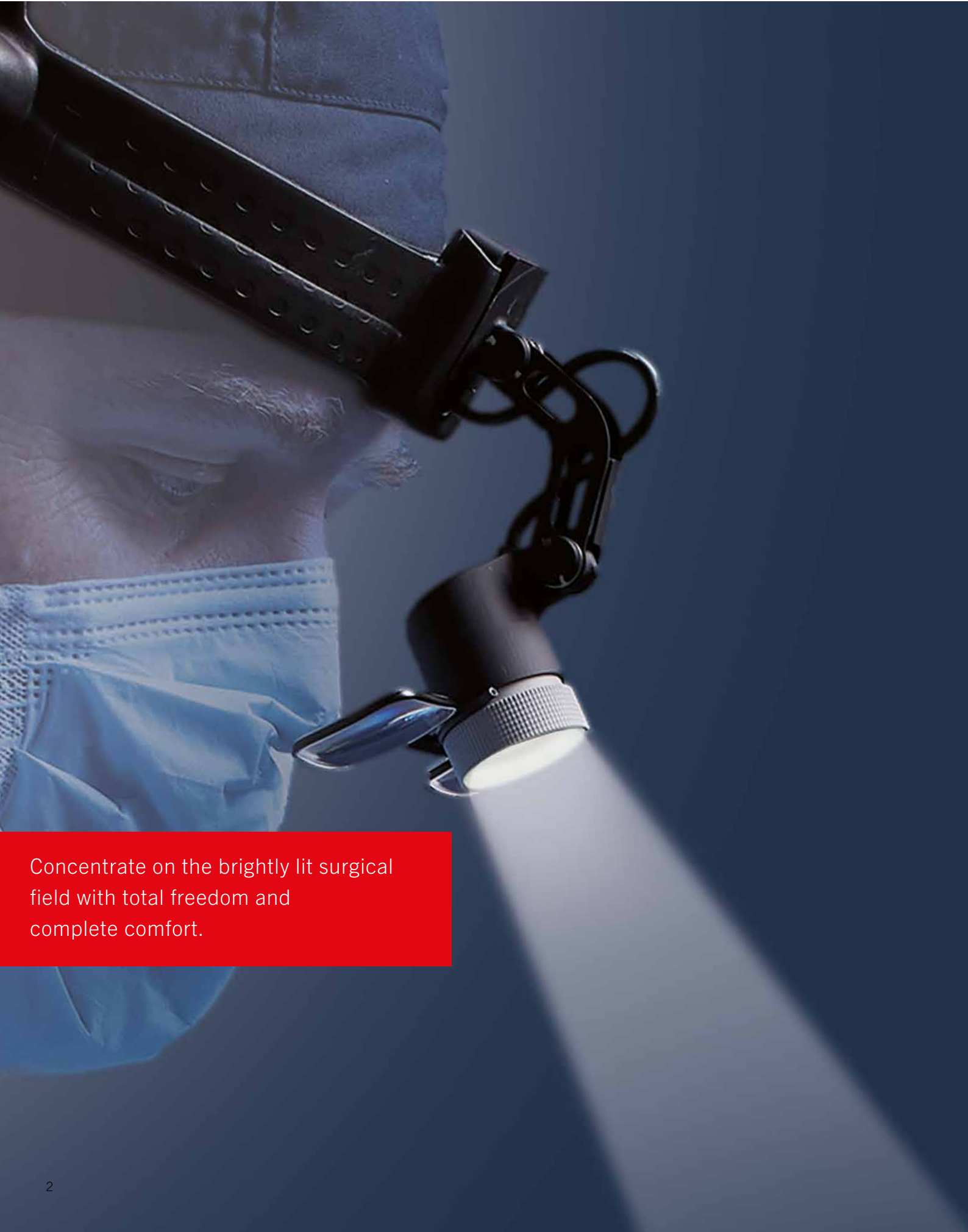

Concentrate on the brightly lit surgical field with total freedom and complete comfort.

## MedLED® Spectra

# Cordless Surgical Headlight

### **Freedom without sacrifice**

The MedLED® Spectra headlight has incorporated the best technology to accomplish the once impossible: a silent, comfortable, cordless headlight that is versatile across many specialties in the clinic or operating room and will enhance your surgical experience in terms of both intensity and spot size.

### **Features luminous intensity:**

- G3 - 100,000 lux
- G6 - 200,000 lux
- G8 - 300,000 lux – now even brighter

### **Battery service life (infinite with 4-bay charger):**

- G3 - 10 hours
- G6 - 8 hours
- G8 - 5 hours

### **Lightweight:**

Headlight weighs 260 g

### **USB charging time:**

4 hours (2 batteries)

### **4-Bay charging time:**

2 hours (4 batteries); charge up to two sets at a time

### **Adjustability:**

Full beam angle adjustability and increased spot size range

### **Spot quality:**

Crisp, circular beam of homogenous light

### **Convertible:**

Interchangeable headband for hard or soft top

### Freedom without sacrifice

If you have a need for a powerful, long lasting, comfortable headlight with no cord, KLS Martin is able to offer you a full range of headlights and accessories.

200 mm  
**G8**

100 mm  
**G3**


90 mm  
**G6**

**maximum spot diameter range**  
Scale 1:1

### Unparalleled brightness

MedLED® Spectra is the light that tailors to you with two headbands and three light options. Choose a light that can power a beam bright enough to see all you need in clinical and OR settings.

### Power and brightness

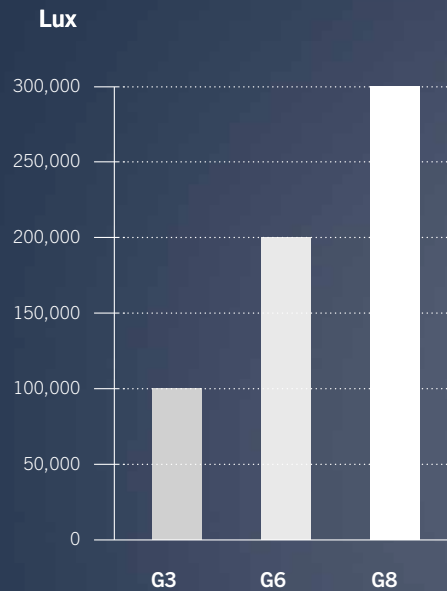

### Quiet-Cool™ Heatsink technology

Cool and quiet  
Proprietary technology is used to cool the LED instantly, avoiding the potential for unnecessary heat buildup. We also designed the housing itself to manage any remaining thermal dissipation.

Fans also demand more power, requiring large batteries attached to the waist of the user. This makes fans a thing of the past.



G3

G6

G8

SPECTRA

Spot adjustability and control



G3



G6



G8

## Convertible Top

Kit comes convertible with both hard top and soft top to choose from for optimal comfort.

Hard top



Soft top



Snap on and off for convenient top conversion




## Go Cordless

Untether from cords and heavy belt battery packs. Dual Everlast™ batteries are integrated into the patented headband and can be changed without losing power.

### Battery life

Hours

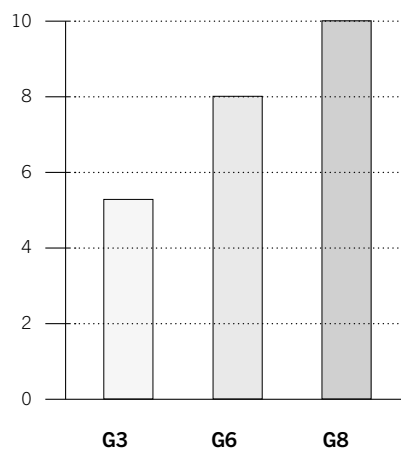

### Extra power

The latest in lithium battery technology



### 4-bay desktop charger

Charge through backup station






### Rapid USB-C charger

Charge through headband



**MedLED®** Spectra  
Headlights  
Ready to use!



	<p><b>MedLED® Spectra G3</b></p> <p>Consisting of:</p> <ul style="list-style-type: none"> <li>■ Spectra G3 headlight</li> <li>■ Carrying case</li> <li>■ Universal power plug</li> <li>■ Padding conversion kit</li> <li>■ 2 rechargeable lithium-polymer batteries</li> </ul>	<p><b>Item No.</b></p> <p>17-510-30-04 Spectra G3 standard kit*</p> <p>17-510-31-04 Spectra G3 hospital kit*</p>
	<p><b>MedLED® Spectra G6</b></p> <p>Consisting of:</p> <ul style="list-style-type: none"> <li>■ Spectra G6 headlight</li> <li>■ Carrying case</li> <li>■ Universal power plug</li> <li>■ Padding conversion kit</li> <li>■ 2 rechargeable lithium-polymer batteries</li> </ul>	<p><b>Item No.</b></p> <p>17-510-60-04 Spectra G6 standard kit*</p> <p>17-510-61-04 Spectra G6 hospital kit*</p>
	<p><b>MedLED® Spectra G8</b></p> <p>Consisting of:</p> <ul style="list-style-type: none"> <li>■ Spectra G8 headlight</li> <li>■ Carrying case</li> <li>■ Universal power plug</li> <li>■ Padding conversion kit</li> <li>■ 2 rechargeable lithium-polymer batteries</li> </ul>	<p><b>Item No.</b></p> <p>17-510-80-04 Spectra G8 standard kit*</p> <p>17-510-81-04 Spectra G8 hospital kit*</p>

\*Hospital Kit (HK) comes with 4 extra batteries and a 4-bay desktop charger (Item No. 17-510-13-04) for backup battery power in a high usage environment.



**MedLED®** Spectra  
Accessories



**Magnifier lenses G3**

	Diopter	Magnification	Increase in eyesight	Optical power	Recommended working distance (focal length)	Item No.
	+ 2	+ 0.40	40%	1.4 x	250 mm	17-500-12-04
	+ 3	+ 0.60	60%	1.6 x	220 mm	17-500-13-04
	+4	+ 0.80	80%	1.8 x	200 mm	17-500-14-04

**Magnifier lenses G6 / G8**

	Diopter	Magnification	Increase in eyesight	Optical power	Recommended working distance (focal length)	Item No.
	+ 2	+ 0.40	40%	1.4 x	250 mm	17-500-22-04
	+ 3	+ 0.60	60%	1.6 x	220 mm	17-500-23-04
	+4	+ 0.80	80%	1.8 x	200 mm	17-500-24-04

## MedLED® Spectra Accessories



17-510-13-04  
4-bay desktop charger



17-510-10-04  
Spare batteries,  
rechargeable (2 pcs.)



17-510-11-04  
USB-C power plug




17-510-06-04  
Spare padding set



17-510-15-04  
Camera adapter



17-500-29-04  
Carrying case, hard

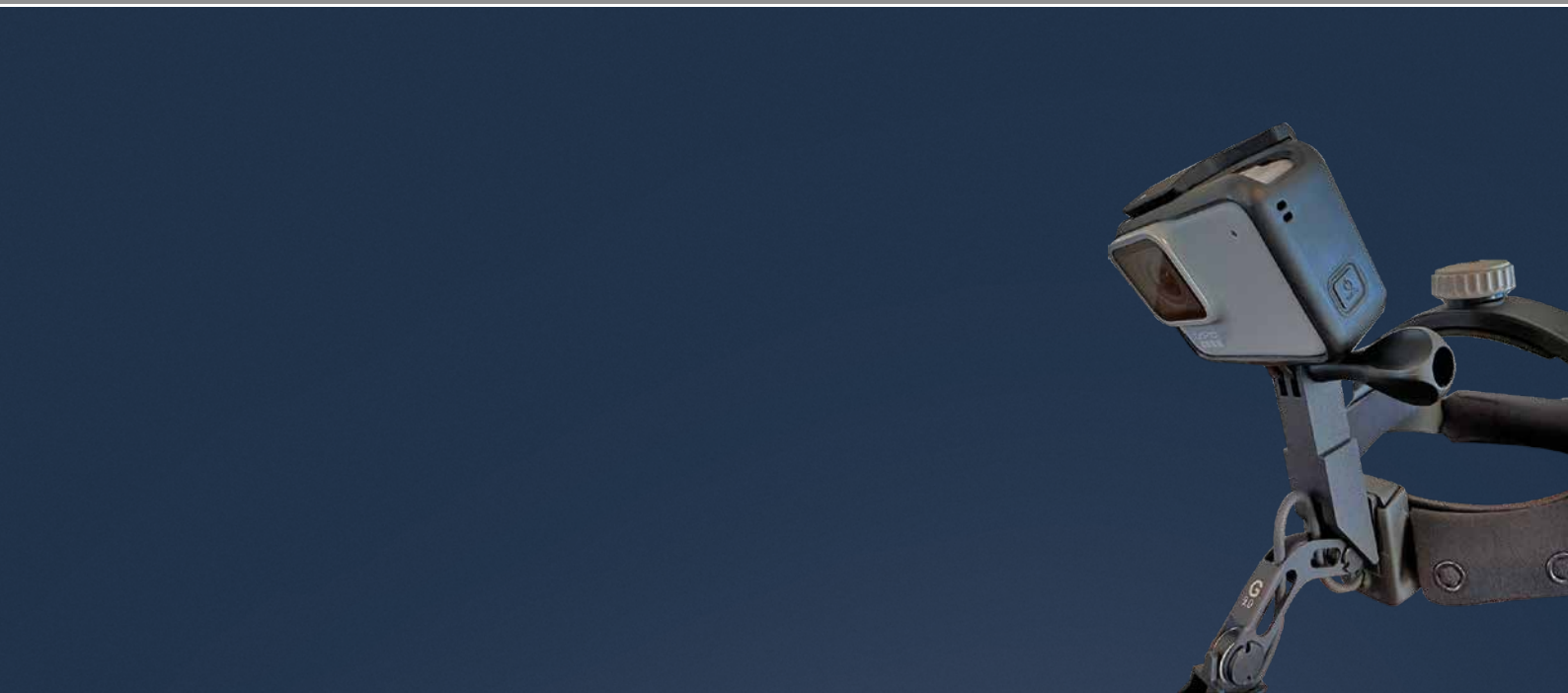

**Bezel**



17-510-32-04  
G3



17-510-62-04  
G6



17-510-82-04  
G8

**Amber filter (for dental)**



17-500-34-04  
G3



17-500-35-04  
G6

17-510-08-04  
Spare head unit

17-510-07-04  
Headband set (HT and ST)



17-510-14-04  
4-bay wall charge station

## KLS Martin Group

### **KLS Martin Australia Pty Ltd.**

Sydney · Australia  
Tel. +61 2 9439 5316  
australia@klsmartin.com

### **Martin Italia S.r.l.**

Milan · Italy  
Tel. +39 039 605 67 31  
italia@klsmartin.com

### **Martin Nederland/Marned B.V.**

Huizen · Netherlands  
Tel. +31 35 523 45 38  
infonl@klsmartin.com

### **KLS Martin UK Ltd.**

Reading · United Kingdom  
Tel. +44 118 467 1500  
uk@klsmartin.com

### **KLS Martin do Brasil Ltda.**

São Paulo · Brazil  
Tel. +55 11 3554 2299  
brazil@klsmartin.com

### **Nippon Martin K.K.**

Tokyo · Japan  
Tel. +81 3 3814 1431  
nippon@klsmartin.com

### **Gebrüder Martin GmbH & Co. KG**

Moscow · Russia  
Tel. +7 499 792 76 19  
russia@klsmartin.com

### **KLS Martin LP**

Jacksonville · Florida, USA  
Tel. +1 904 641 77 46  
usa@klsmartin.com

### **KLS Martin Medical (Shanghai) International Trading Co., Ltd**

Shanghai · China  
Tel. +86 21 5820 6251  
china@klsmartin.com

### **KLS Martin SE Asia Sdn. Bhd.**

Penang · Malaysia  
Tel. +604 506 2380  
malaysia@klsmartin.com

### **KLS Martin Taiwan Ltd.**

Taipei 106 · Taiwan  
Tel. +886 2 2325 3169  
taiwan@klsmartin.com

### **KLS Martin India Pvt Ltd.**

Chennai · India  
Tel. +91 44 66 442 300  
india@klsmartin.com

### **KLS Martin de México S.A. de C.V.**

Mexico City · Mexico  
Tel. +52 55 7572 0944  
mexico@klsmartin.com

### **Gebrüder Martin GmbH & Co. KG**

Dubai · United Arab Emirates  
Tel. +971 4 454 16 55  
middleeast@klsmartin.com

### **Gebrüder Martin GmbH & Co. KG**

#### **A company of the KLS Martin Group**

KLS Martin Platz 1 · 78532 Tuttlingen · Germany  
P.O. Box 60 · 78501 Tuttlingen · Germany  
Tel. +49 7461 706-0 · Fax +49 7461 706-193  
info@klsmartin.com · www.klsmartin.com